Instructions for use

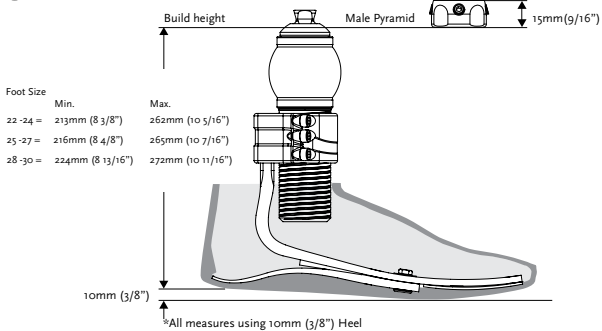
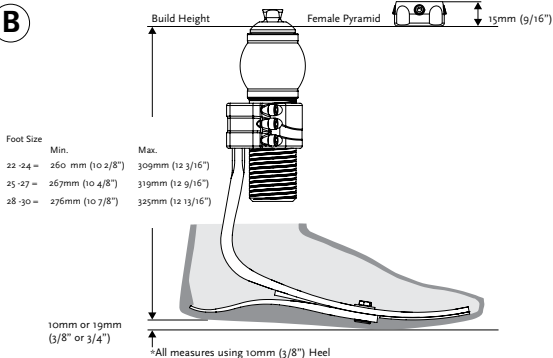
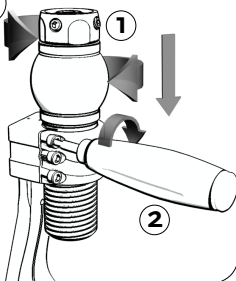
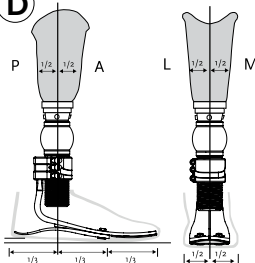
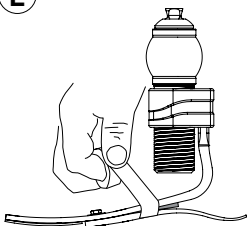
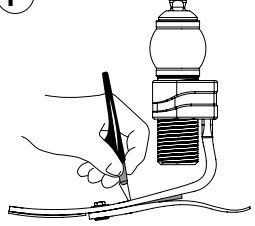
Re-Flex Rotate with EVO™

Life Without Limitations®

.....	3
①EN Instructions for use.....	5
①DE Gebrauchsanweisung.....	7
①FR Guide de fabrication.....	10
①ES Instrucciones para el uso.....	12
①IT Istruzioni per l'uso.....	14
①DA Brugsanvisning.....	16
①SV Bruksanvisning.....	18
①NL Gebruiksaanwijzing.....	20
①PT Instruções de Utilização.....	22
①JA 取扱説明書.....	24
①ZH 中文说明书.....	26

A**B****C****D****E****F**

A complete Re-Flex Rotate foot system comprises the following:

- Re-Flex foot module with EVO
- Maintenance free shock module with rotational compliance
- A cosmetic foot cover with EVO (ordered separately).

Options:

Re-Flex Rotate is available with either a male or female pyramid connector.

Re-Flex foot modules are available either as Standard (**Figure A**) or Tall (**Figure B**).

Warning: Avoid exposure to fire, excessive heat or cold. This may damage the shock or change its function.

Pyramid Adapter

Use Loctite 243 and torque to 15Nm (11ft-lbs).

Flex-Foot Sock with EVO

The toe end of the sock is divided by a thread and fits within the toe split. Upon fitting the sock, hold the top of the toe thread and slide the sock into the toe split. The thread of the sock should rest on top of foot module and not underneath the toe split. The sock should be loose from the heel area.

Shoe Horn

When removing and installing the foot cover use a Flex-Foot shoehorn to avoid damage to foot and/or foot cover.

Alignment Recommendations

To ensure optimum function of the foot and pylon, use the following guidelines:

1. Bench Alignment

To introduce the appropriate socket flexion and heel height:

- Divide the foot length into 3 equal portions (**Figure D**).
- The load line should fall at the junction of the posterior and middle third portions of the foot (**Figure D**).

2. Adjusting the Shock Module (Figure 1)

To adjust the shock module, loosen the three hex screws on the attachment clamp (**Figure C**) and rotate the foot until the correct height and/or rotation is achieved. One full rotation results in 8mm change in height. For finer adjustments, rotate in 90° turns:

Turn:	Build height change
- 90°	2mm lift
- 180°	4mm lift
- 270°	6mm lift
- 360°	8mm lift
+ 90°	2mm lowering
+ 180°	4mm lowering
+ 270°	6mm lowering
+ 360°	8mm lowering

Warning: Ensure there is full engagement of the clamp and shock module. Less than full engagement will void the warranty.

3. **Securing the Shock Module**

To secure the shock module tighten the 3 hex screws (**Figure C,2**). Start with the top hex screw and work down by applying 3Nm with a torque wrench. Perform this sequence twice, since the upper hex screws will loosen in the first run.

Warning: Overtightening the hex screws can strip threads and cause shock module to bind.

4. **Applying a Permanent Shock Module Setting**

To permanently fix the shock module, remove the hex screws one by one; apply Loctite 243 and 3Nm with a torque wrench. When finished, re-apply torque of 3Nm on all hex screws. Failing to apply Loctite 243 to the hex screws may result in the shock coming loose.

Dynamic Alignment

Heel Wedges

Heel wedges influence the heel to toe function. Start altering the heel resistance to improve foot response by adding a heel wedge. The small, medium and large wedges can be interchanged to customize the stiffness and achieve the desired functional characteristics. The wedges can be trimmed using sharp scissors to customize stiffness.

Temporary Wedge Placement (Figure.E)

- Place the wedge in the angle of heel and foot module.
- Secure in position with tape wrapped around the foot module

Permanent Wedge Placement

- Roughen the upper and lower surface of the heel with abrasive paper.
- Apply adhesive on the lower side of the wedge only.
- Locate in the foot/heel junction and position before adhesive sets.
- For split toe feet install the heel wedge then remove a thin slice in the middle by cutting with a sharp knife through the split in the carbon foot module (**Figure F**).

Instant adhesive is necessary to bond the urethane heel wedge. The adhesive cures in 15-20 seconds. For removal the adhesive may be softened by soaking in acetone or cyanoacrylate adhesive remover.

Heel Divider

The heel divider should be placed in the last third of the split heel part. To secure its position a drop of instant adhesive can be used on one side.

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EN – Caution: Össur products and components are designed and tested according to the applicable official standards or an in-house defined standard when no official standard applies. Compatibility and compliance with these standards is achieved only when Össur products are used with other recommended Össur components. If unusual movement or product wear is detected in a structural part of a device at any time, the patient should be instructed to immediately discontinue use of the device and consult his/her clinical specialist. This product has been designed and tested based on single patient usage. This device should NOT be used by multiple patients. If any problems occur with the use of this product, immediately contact your medical professional.

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Össur Iberia S.A.